TRANSCEND[®]

AUTO

Transcend[®] AUTO Quick Guide

Provided for your reference by

cpapXchange.com

ENGLISH

FRANÇAIS

DEUTSCH

ITALIANO

ESPAÑOL

NEDERLANDS

PORTUGUÊS

Note: Transcend AUTO™ is an auto-titrating positive airway pressure device (APAP) for treatment of obstructive sleep apnea. For more detailed device descriptions, cautions and troubleshooting procedures refer to the Transcend AUTO Sleep Apnea Therapy User Manual that is available to download at www.mytranscend.com/resources or by calling Somnetics customer service at 877-621-9626 to request a copy.

WARNINGS

- This device is not intended for water intake and therefore should not be exposed to environmental conditions where the system may get wet.
- This device is not intended for life support.
- The Transcend AUTO must be set up and adjusted by a trained healthcare professional before being used for therapy ramp and pressure.
- The air temperature produced by this device may be as much as 10°F higher than the temperature of the room. Exercise caution if the room temperature is warmer than 90°F (32°C).
- Do not block or otherwise obstruct the exhalation ports of the interface. Follow the instructions for your interface (mask).
- This equipment is not suitable for use with oxygen or in the presence of flammable anesthetic mixtures with air, oxygen or nitrous oxide.
- Use of accessories other than those defined in this manual is potentially unsafe. Use of accessories not recommended herein may result in increased electromagnetic emissions or decreased electromagnetic immunity of the APAP system.
- The APAP system is not defibrillation proof.
- The device should be used only with masks and connectors recommended by Somnetics or recommended by a healthcare professional.

TRANSCEND INTRODUCTION

R_x only

Thank you for choosing Transcend AUTO, the first auto-titrating obstructive sleep apnea therapy device designed to fit your active lifestyle. Weighing less than a pound (426 g) the Transcend AUTO is the smallest, lightest and most innovative APAP system, truly giving you the freedom to sleep anywhere.

Begin by inspecting all items shipped to you. Please report any damaged or missing items to the home healthcare provider where you obtained your device.

Included with the Transcend AUTO:

- Transcend AUTO
- Standard 6 foot (183 cm) hose – not HME compatible (see below)
- Universal hose adaptor
- Multi-plug Universal AC Power Supply Set (PSA2)
- Travel Bag
- CD containing user manuals
- CD containing compliance software
- Quick Guide
- USB cable

What's Not Included (optional accessories sold separately):

- Transcend Base Station (with LCD display or Basic)
- Transcend P₄ Overnight™ and P₈ Multi-night™ Battery Systems
- Transcend Portable Solar Charger
- Transcend Mobile Power Adaptor

OPTIONAL ACCESSORIES

TRANSCEND PORTABLE BATTERY OPTIONS

Transcend P₄ Overnight™ Battery: Weighs 0.6 pounds (272 g). Charge with the Transcend Multi-plug Universal Power Supply, Mobile Power Adaptor, or Portable Solar Charger. 7-9 hours therapy at 14 cmH₂O.

Transcend P₈ Multi-Night™ Battery: Weighs 1.1 pounds (499 g). Charge with the Transcend Multi-plug Universal Power Supply, Mobile Power Adaptor, or Portable Solar Charger. 14-16 hours therapy at 14 cmH₂O.

TRANSCEND DC MOBILE POWER ADAPTOR

The DC Mobile Power Adaptor lets you use Transcend in your vehicle, RV, or boat. The Mobile Power Adaptor may also be used to recharge your Transcend batteries when wall power is unavailable.

TRANSCEND PORTABLE SOLAR BATTERY CHARGER

The Portable Solar Battery Charger uses solar power to recharge your Transcend battery. Solar modules begin charging immediately when placed in the sun, even in low light conditions. Rolled up, it's about the size of a folded newspaper and easily stows for your next outdoor adventure.

OPTIONAL ACCESSORIES

TRANSCEND BASE STATIONS

Transcend Base Station: The Transcend Base Station provides extra stability as a bedside anchor for the Transcend CPAP. While not required for home use of the Transcend CPAP, in some cases use of the Base Station may be preferred.

Transcend LCD Programmable Base Station:

In addition to providing enhanced bedside stability, the Transcend LCD Programmable Base Station allows users to view basic therapy compliance information. The Programmable Base Station may also be used by the clinician to change therapy pressure and ramp settings when used in the Clinician Mode (Clinician Mode access is not available to patient users).

ASSEMBLING YOUR TRANSCEND AUTO

Transcend **AUTO** Sleep Apnea Therapy System

To completely assemble the Transcend Sleep Apnea Therapy System:

1. Attach the Universal Hose Adaptor to the 6-foot hose;
2. Attach the Universal Hose Adaptor and 6-foot hose to the APAP;
3. Connect the interface (mask) to the patient (open) end of the 6-foot hose;
4. Plug the power supply barrel connector into the APAP with the cord pointed up;
5. Connect the power supply to a power source.


POWERING YOUR TRANSCEND

Your Transcend comes with a Multi-plug Universal AC Power Supply set that is used with a wall power outlet. The power supply can be used worldwide with the included exchangeable plug adaptors. Do not use a power converter with your Transcend as it may damage the device.

To use the Multi-plug Universal AC Power Supply:

1. The Multi-plug Universal AC Power Supply comes with several plug blades for global use (Image 1). If the correct plug blade is not already attached to the power supply, depress the button on the plug blade and turn the plug blade counterclockwise until it releases from the power supply.
2. To attach the desired plug blade, line up the protruding tabs on the inside of the plug blade with the slots on the power supply and gently push the tabs into the slots. Turn clockwise until the blade clicks into place.
3. Insert the barrel connector of the power supply into the power jack on the back side of the device (image 2). Be sure the cord is facing upward when connected to the device.
4. Connect the wall mount power supply plug to a wall power outlet.
5. The control panel of the Transcend AUTO has green and yellow LED indicator lights. When power is supplied to the APAP, the yellow LED light will flash once followed by three quick flashes of the green LED light. This sequence indicates power is being supplied to the device and that it has successfully entered Standby mode.


Your Transcend AUTO can also be powered by the optional Transcend Mobile Power Adaptor or Transcend Battery Systems.

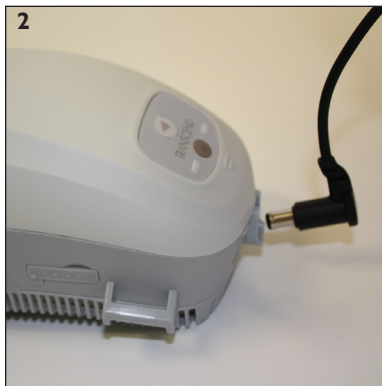STARTING AND ENDING THERAPY

Starting Therapy

1. The Transcend control panel consists of a Power button (bottom), Ramp button (top), and two LED lights. (Refer to image)
2. Ensure that the device is connected to a power source and is in Standby mode. (See Powering Your Transcend)

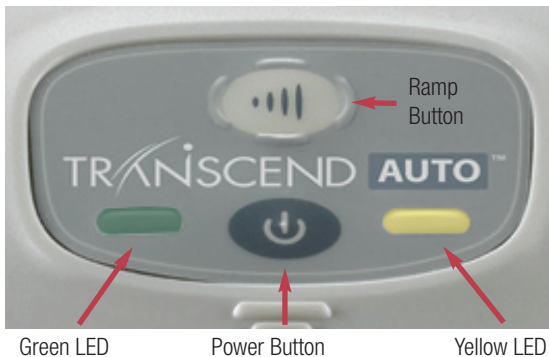3. Put on your mask and make sure that proper fit is achieved.
4. Press the lower, gray Power button () to turn on Transcend and enter Therapy mode. Air flow will begin as the blower delivers your prescribed therapy pressure.

Using the Ramp Function

Ramp settings will be programmed by your healthcare provider as appropriate.

1. Initiate ramp by pressing the Ramp (top) button after the Power button has been pressed.

Note: The Ramp button on the Transcend will not work unless your healthcare provider has preset the ramp function.

2. Pressure from the device will be lowered to a preset starting pressure and will gradually increase over a preset time interval until the device reaches your prescribed starting therapy pressure.
3. To accelerate the ramp time, hold down the Ramp button. The device will quickly increase pressure until reaching your prescribed therapy pressure.

Ending Therapy

To end therapy, press the Power button to return the device to Standby mode.

CARING FOR YOUR TRANSCEND AUTO AND COMPONENTS

Drying Mode

Drying mode should be used **daily** to dry the inside of your Transcend after use.

Caution: Do not attempt to use Transcend while in Drying mode. The device will not provide prescribed therapy pressure while in Drying mode.

1. Set the device on a dry, stable surface. A cloth surface is not recommended as it may get soiled.
2. Make sure the device is in Standby mode.
3. To enable Drying mode, press and hold the Ramp button and then press the Power button.
4. Low pressure airflow will begin and last for 30 minutes until the drying sequence is completed. The device will then return to Standby mode. To end drying mode before the 30 minute sequence is completed, press the Power button. The device will return to Standby mode.

Caution: Drying mode will dry and clean the inside of your Transcend. Do not attempt to use any other method to clean the inside of the Transcend device.

Cleaning the Exterior

Clean the exterior of your Transcend device as needed. To clean, wipe the surface with a lint-free cloth or paper towel that has been dampened in water containing mild dishwashing detergent. Following cleaning, wipe again with a cloth or paper towel moistened with clear water. Allow to dry.

Caution: Never submerge Transcend in water. Do not allow soap or water to enter the device.

Cleaning the Transcend Components

The interface (mask) should be washed weekly. The air supply tube (hose) should be washed daily.

1. Completely disconnect the components.
2. Hand wash in warm water with a mild dishwashing detergent.
3. Rinse thoroughly.
4. Allow to air dry. Components should be completely dry before reconnecting.
5. Reconnect the components and re-attach to Transcend.



Caution: Do not put any Transcend components in the dishwasher or laundry machine.

CLEANING AND REPLACING THE AIR INLET FILTER

The Transcend AUTO Air Inlet Filter should be replaced every six months. Contact your homecare dealer to purchase a replacement Air Inlet Filter.

1. Remove the filter plate by gently pulling back the connector and lifting the plate away from the APAP (Image 1).
2. Remove the filter from the filter frame and discard.
3. Insert a new filter into the filter frame.
4. To replace the new filter plate, lay it on the bottom of the APAP. Snap the front clips into place first, then snap the connector into place (Image 2).

The Air Inlet Filter should be removed and cleaned weekly. To clean the Filter, remove from the device, wash thoroughly in water containing mild dishwasher detergent, rinse in clear water, and air dry. The Filter should be completely dry before returning it to the device. Replace the Filter as described above.

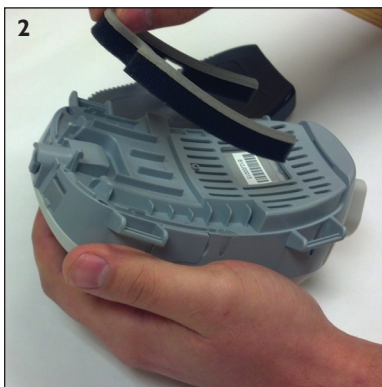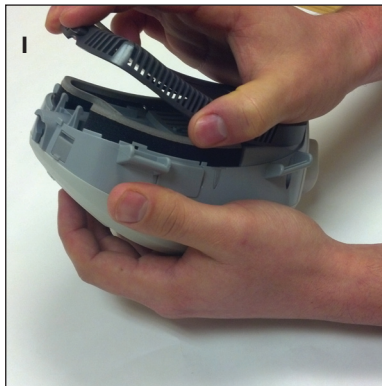Note: Earlier versions of Transcend may require that the complete filter frame be cleaned and replaced. Refer to the User Manual supplied with the device.

SOMNETICS INTERNATIONAL, INC.

33 5TH AVE NW, SUITE 500
NEW BRIGHTON, MN 55112
PHONE: 651.621.1800
TOLL FREE: 877.621.9626
FAX: 651.204.0064
www.mytranscend.com

AUTHORIZED REPRESENTATIVE (MDD 93/42/EEC)

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